

SpeediCath® catheter

Single use hydrophilic coated catheter for intermittent catheterization
– comes ready to use in a sterile saline solution

- Ready for instant use
- Extremely smooth surface
- Hygienic to handle
- New advance in coated catheters



SpeediCath[®] catheter

The SpeediCath family is an innovative range of single-use hydrophilic coated catheters for intermittent catheterization


Ready for instant use

SpeediCath is a PVC-free hydrophilic coated catheter. It comes packaged in a saline solution, making it ready for immediate use when you open the packaging. There is no need for water, nor must you wait the usual 30 seconds for the coating to be activated.



Extremely smooth surface

The hydrophilic coating covers not only the outside surface but even the inside of the catheter eyelets. An exclusive new production process gives the coating a unique ability to bind and retain water, making it exceptionally smooth and even.

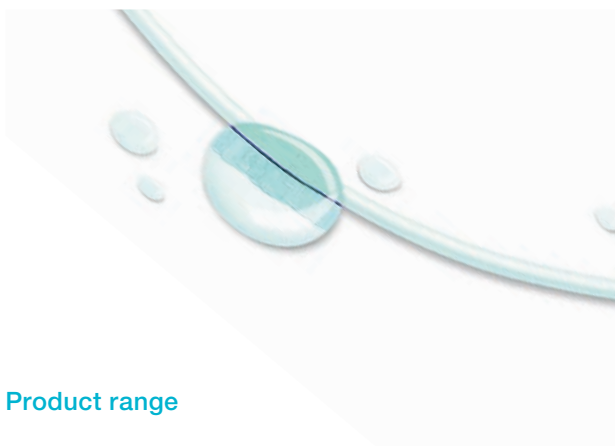

Product range

SpeediCath Catheter



Size	Men Nelaton	Men Tiemann	Female	Boys	Paediatric
06			28506	28606	28706
08	28408		28508	28608	28708
10	28410	28490	28510	28610	28710
12	28412	28492	28512	28612	
14	28414	28494	28514		
16	28416		28516		
18	28418				

The SpeediCath catheter is for single use only. Units per box: 30

Ostomy Care
Urology & Continence Care
Wound & Skin Care

The smoothness dramatically reduces friction between the catheter and the urethral mucosa and minimises any risk of long term damage.



Hygienic handling

The packaging which can be opened at both ends, can be held at one end (acting like a glove) while the catheter is inserted from the other. This makes it more hygienic for users as well as for care givers.

Discreet solution

The SpeediCath catheter is foldable and can be carried discreetly in a pocket or small handbag. However, it should not be kept folded for more than 3 hours.

A complete system

The range also includes Conveen urisheaths and urine bags. Together with EasiCath and SpeediCath catheters, this is a complete system of high-quality appliances for people with bladder emptying problems.

Composition

Thermoplastic polyurethane, polyvinylpyrrolidene, water.



Coloplast develops products and services that make life easier for people with very personal and private medical conditions. Working closely with the people who use our products, we create solutions that are sensitive to their special needs. We call this intimate healthcare. Our business includes ostomy care, urology and continence care and wound and skin care. We operate globally and employ more than 7,000 people.

The Coloplast logo is a registered trademark of Coloplast A/S. © [2011-08.]
All rights reserved Coloplast A/S, 3050 Humlebæk, Denmark.

Coloplast A/S
Holtedam 1
3050 Humlebæk
Denmark

www.coloplast.com